

POCT Oral Fluid Testing

Choose one assigned POCT device • follow the panel printed on that device • use Quantisal only if laboratory confirmation is required

BEST PRACTICE • DAY-OF-COLLECTION CHECKLIST

After countersigning the consent form, follow the interactive checklist step by step. It is the best way to ensure the correct process and protocols are completed.

FIELD USE

One donor, one POCT device. Confirm the assigned device before you begin. The BTNX Rapid Response Saliva Test Cup and the BTNX Oral Fluid Lollipop are **alternative screening options**; do not use both on the same donor.

OPTION 1 • SALIVA TEST CUP

When: The collection is assigned a BTNX Rapid Response Saliva Test Cup.

Action: Use the cup pathway on pages 2–3. Follow the drug panel and timing printed on the supplied cup and its insert.

OPTION 2 • ORAL-FLUID LOLLIPOP

When: The collection is assigned a BTNX oral-fluid lollipop.

Action: Use the lollipop pathway on pages 4–5. Follow the exact drug codes, cut-offs, and result window printed on the supplied lollipop and its insert.

Use the Saliva Test Cup

Follow this page only when the BTNX Rapid Response Saliva Test Cup is the assigned POCT device

PRODUCT AT A GLANCE

BTNX Rapid Response Saliva Test Cup

Use the sealed cup supplied for the collection. Follow the panel labels and product insert for the specific cup configuration assigned to this donor.

Photo: BTNX Inc.



PRE-COLLECTION CONTROL

10 min

Before either POCT device is opened, nothing may enter the mouth: no food, drink, gum, tobacco or nicotine, lozenge, spray, mouthwash, or other material.

1

Check and open the assigned cup

Confirm the cup package is sealed and in date. Record the lot number and expiry in the required collection record before opening. Open the cup in the donor's presence.

2

Collect active oral fluid

Keep the collection in view. Have the donor actively swab the cheeks and gumline. Do not use passive placement or leave the donor alone during collection.

Read the Saliva Test Cup

Complete saturation • seal and label • then use the timing printed on the assigned cup

3

Wait for full saturation

Do not remove the swab until the saturation indicator is solid and uniform red. A partial or faint indicator is not sufficient.

4

Secure and start the timer

Have the donor place the swab in the collection tube. Seal and label it immediately, record the start time and collector ID, then set the timer. Read the result only in the window stated on the supplied cup and insert.

Negative



See: For a drug strip: control line present and any test line present—even a faint line.

Do: Document the negative screen in the ScreenClear record. Do not collect Quantisal unless laboratory confirmation is required.

Non-Negative



See: For a drug strip: control line present and no test line.

Do: Presumptive only. Record the affected drug code(s); do not call the result positive. If laboratory confirmation is required, collect Quantisal during the same encounter and begin the Dynacare workflow.

Use the Oral-Fluid Lollipop

Follow this page only when a BTNX oral-fluid lollipop is the assigned POCT device

PRODUCT AT A GLANCE

BTNX Oral Fluid Lollipop

Use the lollipop supplied for the collection. It is a qualitative POCT screen only—not a final laboratory result. Follow its label and insert for the assigned panel.

Photo: BTNX Inc.



CHECK THE ASSIGNED PANEL

FIELD USE

Before opening the lollipop, verify the exact drug codes and screen cut-offs printed on the supplied device and insert. Different lollipop configurations can have different drug panels. Record the assigned panel, lot number, and expiry in the required collection record.



DO NOT SUBSTITUTE A PANEL

Use only the lollipop configuration assigned in the interactive Day-of-Collection Checklist. Do not choose, add, or remove a drug panel independently.

Collect & Read the Lollipop

Keep the device level • wait for control lines • use the result window printed on the assigned lollipop

5

Protect the mouth-clear period

For at least **10 minutes**, nothing may enter the mouth. If anything enters the mouth, restart the full 10-minute period and document the restart.

6

Collect oral fluid

Open the pouch only when ready. Remove the cap and place the collector in the donor's mouth. Keep the device level until saliva flows across the test strips and the control lines appear.

7

Cap and time the lollipop

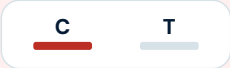
Replace the cap, place the device on a level surface, and start the timer immediately. Do not tilt, open, or handle the result window while waiting.

DRUG RESULT WINDOW

Use the label

Follow the read window printed on the assigned lollipop and in its insert. Do not interpret a result outside that window.

Non-Negative

C

T

See: For a drug strip: control line present and no test line for a code on the assigned lollipop panel.

Do: Presumptive only. Record the exact affected lollipop code(s). If laboratory confirmation is required, collect Quantisal during the same encounter and begin the Dynacare workflow.

If Confirmation Is Required

Quantisal is used only after a presumptive non-negative POCT result when laboratory confirmation is required

FIELD USE

Decision point: The used POCT cup or lollipop is not the Dynacare confirmation specimen. Collect Quantisal only when laboratory confirmation is required after a presumptive non-negative POCT result.

MATCH EVERY FLAGGED CODE

1

Record the exact code

Record each non-negative code printed on the assigned cup or lollipop panel.

2

Use the matching confirmation order

Select only the Dynacare laboratory confirmation target(s) that match the flagged code(s) for that assigned device configuration.

3

Do not guess or over-order

Do not add unrelated drug classes, substitute a panel, or make an independent final result determination.

!

PANEL-SPECIFIC CONFIRMATION

The assigned device label and the interactive Day-of-Collection Checklist identify the panel and matching confirmation pathway. If the code or order is unclear, contact ScreenClear before collecting Quantisal.

Collect Quantisal & Send

Use Quantisal only for required laboratory confirmation after a presumptive POCT non-negative

PRODUCT AT A GLANCE

Quantisal Oral Fluid Collection Kit

Use the Dynacare-supplied Quantisal kit only when laboratory confirmation is required. It is the separate collection, preservation, and transport device—not the used POCT cup or lollipop.

Photo: Abbott Toxicology

**8**

Collect to the blue indicator

Check the package and transport tube(s) are intact and in date. Have the donor pool saliva, then place the collector under the tongue with head down. Do not allow chewing, sucking, talking, or removal during collection. Wait until the required indicator(s) are completely blue.

9

Secure, label and pack

Place the collector in its matching transport tube and snap the cap fully closed. For a dual-pad kit, follow the designed break point and matching A/B tubes. Complete the Dynacare CCF/eCCF with the flagged confirmation target(s), then use the supplied specimen bag and document pouch.

!

NO BLUE BY 10 MINUTES

Discard the incomplete Quantisal collector and begin again with a new kit only after saliva has accumulated. Do not use a partial specimen.

Finish the Record

Protect POCT evidence, Quantisal custody, and the required documentation

10 • COMPLETE THE DYNACARE HANDOFF

When: A required Quantisal confirmation specimen and the CCF/eCCF are complete.

Action: Use the approved Purolator process for prompt shipment. Photograph CCF Copy 2 and Copy 4; images upload automatically and are sent to the required destination. Place CCF Copy 3 and required originals in the ScreenClear document folder.



INVALID POCT

No control line means do not interpret the drug strip. Repeat with a new device of the same assigned type. If a valid screen cannot be obtained, document the facts and contact ScreenClear before closing the record.



NO EMPLOYER COMMUNICATION

Do not report any POCT result to the employer or DER. A non-negative is presumptive, not a final laboratory result. Preserve the device and let ScreenClear manage communication.



The interactive Day-of-Collection Checklist is complete.



Exactly one assigned POCT device was used: Saliva Test Cup or Oral-Fluid Lollipop.



If confirmation was required, Quantisal, Dynacare custody, automatic image upload, and document-folder requirements are complete.



STOP • SECURE • DOCUMENT • CONTACT SCREENCLEAR

1-844-925-3274 | info@screenclear.ca
Never contact the employer or DER directly.